

SYSTEM TWO
HALOGEN SERIES

Bovie[®]

SYSTEM TWO
HALOGEN Series

MULTI-FUNCTIONAL SURGERY LIGHTING



CONVENIENT

The SYSTEM TWO Halogen Series by Medical Illumination International offers multi-purpose configurations for your convenience

BRIGHT

The System Two Halogen Series produces up to 102,000 Lux* (9,500 foot candles*) at 1 meter

WHITE

The Halogen series has 4,000° Kelvin color temperature which provides excellent tissue rendition

CLEAR

Multi-faceted reflector design increases shadow reduction for enhanced visibility during surgery

VERSATILE

Two light-head diameters options (a 16" reflector and/or a 20" reflector), four-Stage dimming, a variety of mounting styles and a multitude of optional accessories offer the ultimate degree of office and surgery room capability

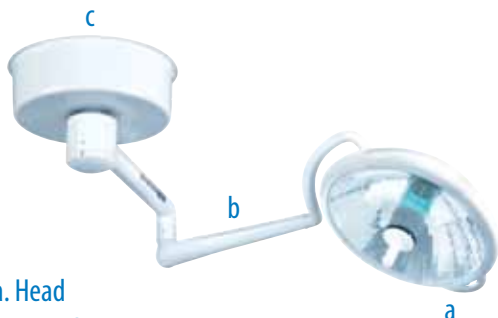
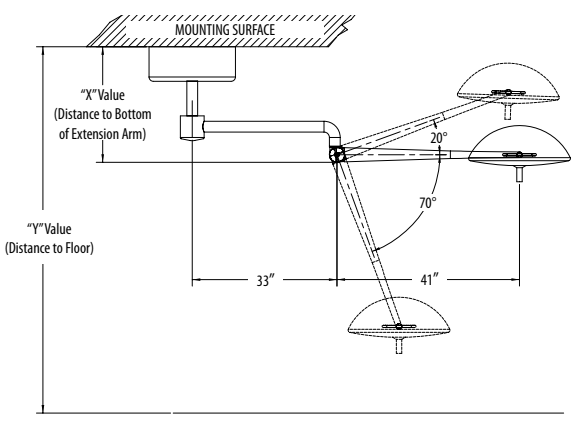
RELIABILITY

The System Two Halogen Series is backed by a three (3) year manufacturer's limited warranty






*Measurement based on S2S-20 light-head.

SURGICAL LIGHTING

SYSTEM TWO
HALOGEN SERIES
INCLUDES

Single Ceiling Mount Model	Rotation and Arm Movement
 <p>a. Head b. Articulating Arm c. Single Ceiling Mount</p>	
S2S-16 (16" Reflector) & S2S-20 (20" Reflector)	

OTHER SYSTEM TWO
HALOGEN MODELS
AVAILABLE

Wall Mount Model	Duo Ceiling Model	Portable Floor Model
 <p>S2WM-16 (16") or S2WM-20 (20")</p>	 <p>S2D-20 (20") & S2D-20 (20")</p>	 <p>S2FM-16 (16") & S2FM-20 (20")</p>
Wall Mount Monitor	Trio Ceiling Mount w/Video Camera	
 <p>S2WM-MA</p>	 <p>S2T-20VC (20") & S2T-16VC (16") *Controller included with Video Camera</p>	

SYSTEM TWO
HALOGEN SERIES
TECHNICAL SPECIFICATIONS

Dimensions: S2S-20[†] Ceiling casting 17.0" (431.8 mm) Dia x 4.7" (119.4 mm) deep Ceiling rod 2.5" (63.5 mm) Dia x 9.0"-43.0" (228.6 - 1092.2 mm) long (depending on ceiling height) Arm Assembly (to center of light-head) 74" (1,880 mm) Light head assembly 23.0" (584.2 mm) diameter x 9.5" (241.3 mm) deep		Rotations*: Ceiling mount/extension arm interface Continuous Articulating/extension arm interface Continuous Articulating arm vertical movement +20, -70 Degrees Articulating arm/Yoke interface Continuous Yoke/lamp head interface +/- 150 Degrees	
Weight*: Solo ceiling mount assembly Approximately 73.5 lbs (33.3 Kg) Duo ceiling mount assembly Approximately 104.5 lbs (47.4 Kg) Trio ceiling mount assembly Approximately 155.5 lbs (70.5 Kg) Light head assembly Approximately 17.5 lbs (7.9 Kg)		Data*: Light Source Halogen Lamp Bulb Life in Hours 1,000 hours (average) Color Temperature 4000 Kelvin Lux 102,000 Lux* (9,500 foot candles*) Input Voltage 100 - 240 VAC 50/60 Hz	
[†] For the S2S-16 model technical specifications, contact your local Bovie Sales Rep.		*Measurement based on S2S-20 light-head	

 **Made in the USA. The SYSTEM TWO Halogen Series is currently only available for sale in the USA.**

